

SPECIFICATION & DATA SHEET


N.L.L.V. (N.I.S.T. LOCKABLE LINE VALVE)

The MIPS N.L.L.V. is designed for economy, maximum design life and reliability whilst requiring minimum maintenance. They are available as a standard product with the following design features:

- Comply fully with the requirements of HTM 02-01
- All valves fitted with BS 5682 N.I.S.T. connectors
- Sealed shut-off N.I.S.T. valves
- All valves fitted with flat-faced 'o' ring unions
- Standard pipework sizes 15-159mm
- Each valve supplied with blanking facility
- Every N.L.L.V. pressure tested prior to despatch
- Every valve tested for tightness (no pressure pass when closed)

In addition to all of this, MIPS can **custom build** N.L.L.V. to your specific requirements including;

- Multi-Valve Manifolds complete with fixing frame
- Padlocks

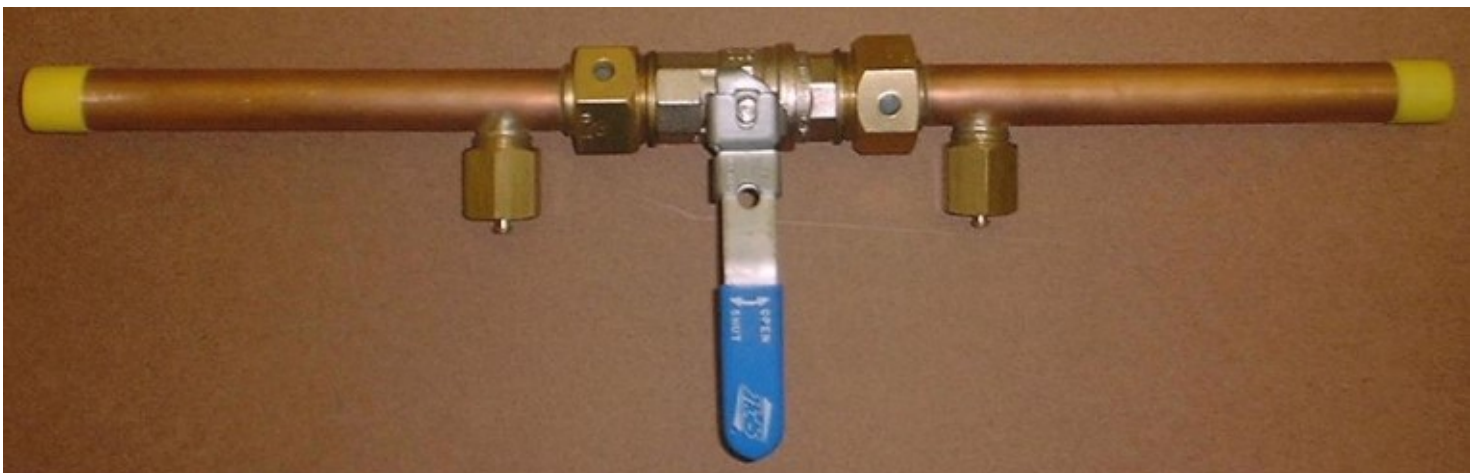

The right is reserved to change the specification of the items described herein at any time.